

Correlative Microscopy - Algorithm overview

Martin Schorb, martin.schorb@embl.de

April 25, 2013

This document can be found in your scripts directory, subfolder doc/

corresponds to commit 5550fbc committed at 2013-04-25 14:43:18 +0200
of martin_correlate.m at <http://git.embl.de/cgit/schorb/corr/>

Contents

1	Introduction	2
2	Installation and requirements	2
2.1	System and software requirements	2
2.2	Obtaining the newest version of the correlation algorithms	2
2.3	Initialization - the init script	4
3	Correlation from LightMicroscopy to LowMag EM image	6
3.1	Executing the script	6
3.2	Output and generated files	7
3.3	User interaction and key procedures	8
3.3.1	Initialization dialog	8
3.3.2	Contrast adjustment - optional	8
3.3.3	Rotation of the fluorescence fiducial image	8
3.3.4	Fiducial selection	9
3.3.5	Fiducial sub-pixel fitting and display	9
3.3.6	Picking the fluorescent spot(s) of interest	11
3.3.7	Sub-pixel fitting of the fluorescent spot(s) of interest	11
3.3.8	Correction of stage drift in between imaging of the different fluorescence channels	13
3.3.9	Calculation of the optimal transform and predicting fiducial and fluorescent spot EM-coordinates.	13
4	Correlation form low magnification tomogram to high magnification EM image	16
4.1	Executing the script	16
4.2	Output and generated files	16
4.3	User interaction and key procedures	17
4.3.1	Fiducial selection	17
4.3.2	Correlation	17
5	Useful hints and tricks	18
5.1	Overlay of correlated images	18

1 Introduction

This document describes the usage of the algorithms designed for correlation of fluorescence microscopy images to their corresponding EM image.

The correlation procedure, as described in Kukulski et al., JCB (2011) and in more detail in Kukulski et al., MethCellBiol (2012) consists of two major coordinate transformations. The first one uses fluorescent microspheres to calculate the mapping of a single point fluorescent signal onto an electron microscopy image of low resolution (4 – 10k \times). The obtained coordinates can then be transformed further using a different fiducial system (in this case the gold beads used for tomogram reconstruction) to a high magnification image.

This document provides a step by step manual in how to use the algorithms, their parameters and outputs and shows the corresponding code snippets to give an idea of the points at which certain things happen while running the script.

The pdf-file of this documentation corresponding to the version of the scripts you are using can be found in your scripts directory, subfolder doc/.

2 Installation and requirements

2.1 System and software requirements

A MATLAB[®] installation version 7.4.0 onwards including the Image Processing Toolbox is required for successfully running the scripts. The Optimization Toolbox is required for sub-pixel Gaussian fitting of the fluorescent spots but the correlation procedure can be run without it.

The scripts should run independent from the type of operating system used. However its function was so far only successfully tested on Linux and Mac OS X 10.4 environments.

In order to access the newest version of the scripts, a Git¹ client is needed. It is usually included in most Linux distributions and easy to obtain for other operating systems. In case you don't have access to a software capable of checking out subversion repositories, you can simply download a snapshot of the scripts directory at <http://git.embl.de/gitweb/?p=schorb/corr.git;a=snapshot;h=master> or all its files individually from the webserver (<http://git.embl.de/cgit/schorb/corr/tree/>, not recommended).

2.2 Obtaining the newest version of the correlation algorithms

The main scripts as well as all supporting underlying algorithms and this documentation are stored in a central Git repository which is accessible from within the EMBL network.

First installation

To get the most recent version of algorithms and documentation files for the first time, you need to copy the code from the server into the directory where you want to store these files. Open a Unix Terminal shell and change to the directory in which you want the scripts to end up. To get the scripts type:

```
git clone git://git.embl.de/schorb/corr.git XXX
```

Where XXX is an optional directory name where the scripts will be copied into. If no directory name is specified a folder `corr/` containing all files will be generated. Accept the encryption certificate.

You should see a result similar to this:

¹See <http://git.embl.de/>, [http://en.wikipedia.org/wiki/Git_\(software\)](http://en.wikipedia.org/wiki/Git_(software))

```
% git clone git://git.embl.de/schorb/corr.git XXX
Cloning into 'XXX'...
remote: Counting objects: 61, done.
remote: Compressing objects: 100% (60/60), done.
remote: Total 61 (delta 7), reused 0 (delta 0)
Receiving objects: 100% (61/61), 4.79 MiB | 9.05 MiB/s, done.
Resolving deltas: 100% (7/7), done.
```

Now you have a local working copy of the most recent versions of the correlation files.
[You can also now find the most up-to-date copy of this manual in the doc/ subfolder.](#)

Add this directory including all subdirectories to the MATLAB path by adding the following line to MATLAB's `startup.m` script, that is usually located in a MATLAB related folder within your home directory (`~/.matlab/startup.m`) when working in a Linux environment.

```
addpath (genpath('/path/to/your/directory'), '-begin')
```

On Macintosh or Windows systems you have to add the path by using MATLAB's preferences. Make sure that also the subdirectories are added.

Updating your existing scripts

To have your scripts and documentation files always up to date according to the current revision, update them from the repository by simply executing

```
git pull
```

in the directory in which you put the scripts.

A successful update of the scripts will look like this:

```
remote: Counting objects: 9, done.
...
Updating f98710a..d4e01df
...
doc/Correlate.tex |      8 ++++++--
2 files changed, 7 insertions(+), 1 deletion(-)
```

If you get an error stating: `fatal: Not a git repository ...`
 Make sure your active terminal is in the correct directory.

If you receive another error message, this most likely means that you've edited some of the scripts yourself and thus there's a conflict resulting, because the file you want to update from the server is different from what the program expects. The easiest solution to this is save your initialization script (see 2.3) to another directory and delete the entire script directory (including all hidden files). Then perform a clean new installation of the scripts as described above.

Editing the scripts

You are able manipulate and edit your local working copy of the scripts (activate or remove sub-pixel fitting, shift correction etc.). However these changes will be reverted while updating to a newer revision. **Therefore please be really cautious while doing this! Use the initialization script instead to personalize the correlation procedure or let me know if there's substantial changes that you need.**

Also note that there will be some hidden files and directories (`.git` etc.) written while checking out. These contain important information needed by Git and should not be modified or deleted.

2.3 Initialization - the init script

You can adjust key parameters of the correlation scripts by modifying the initialization script as shown below. When checking out the repository a file called `corr_init_orig.m` will appear. You can simply rename this file to `corr_init.m` and change the parameters and paths according to your needs. This file will then stay in your local scripts directory and will not be overwritten while updating the scripts with the newest version from the repository.

It might occur that in future revisions additional parameters will be added to this script, so in case you run into an error message stating this, just update your local `corr_init.m` with the changes you find in the downloaded and most up-to-date `corr_init_orig.m`.

In the init script you can specify the following settings:

- directories or directory roots where you typically store your correlations and thus the fiducial, shift and highmag coordinate files.
- whether or not the fluorescence image is flipped compared to the EM image. (This depends on the orientation of the grid in the EM)
- manual contrast settings for the FM images
- skipping of shift correction
- size of the prediction circle
- production of overlay fluorescence images
- correlate multiple or single spots
- sub-pixel localization of fluorescent spots (fiducials and/or POI, this requires the Optimization Toolbox to be installed)
- interactive mode of fitting

```

1 % corr_init()
2 % -----
3 % version MartinSchorb 130312
4 % Copyright EMBL 2013, All rights reserved
5 %
6 % Initializes pathnames and parameters for LM/EM Correlation script
7 % martin_correlate.m
8
9
10 % ----- replace > pwd < by the directory of choice in parentheses -----
11
12 % Location to search for fiducial coordinate file - *.pickspots.mat
13 loc_pickspots = pwd; % '/directory/to/correlation';
14
15 % Location to search for shift correction coordinate file - *.shiftcoos.mat
16
17 loc_shiftcoos = pwd; % '/directory/to/correlation';
18
19 % Location to search for Highmag Fiducial coordinate file - *.lmhmcoos.mat
20
21 loc_hmcoos = pwd; % '/directory/to/HMcorrelation/';
22
23 % Flip fluorescence images (different grid orientation in LM and EM)
24 % flip = 1 if images should be flipped // flip = 0 if images are in same orientation
25
26 flip = 0;
27
28 % Adjust the contrast of display of the fluorescence images (blue, green, red) (0 - auto; 1 - open
29 % adjustment window)
30
31 contr_fid = 1;
32 contr_poi = 0;
33 % obsolete contr_other = 0;
34
35 % Skip the shift adjustment between channels (active if 1)
36
37 shift_skip = 0;
38
39 % Size of prediction circle in nanometers
40
41 accuracy=50;
42
43 % write overlay images for high-mag correlation (files will be written if 1)
44
45 hm_overlays = 0;
46
47 % multiple spots of interest?
48
49 multispot = 0;
50
51 % subpixel localization of fluorophores enable/disable (Optimization Toolbox required!!)
52 % options:
53 % 0 - no fitting
54 % 1 - fit signal of interest only
55 % 2 - fit fiducials only
56 % 3 - fit both fiducials and signal of interest
57
58 gaussloc = 3;
59
60 % interactive mode for fitting (0 - inactive, 1 - active, 2 - always active)
61
62 fit_interactive = 2;
63
64 % -----
65
66 % other parameters that can be adjusted, if unsure leave them as they are
67
68 trafo = 'linear_conformal';
69 % trafo = 'affine';
70 % trafo = 'projective';
71
72 fmbboxsize=11; % must be odd number - size of box for bead-image subpixel fitting
73 imbboxsize=11; % must be odd number - size of box for fluo-image subpixel fitting
74
75 fmfilt = 45;
76
77 pixelsize_lm = 5.068; % pixel size lowmag tomogram in nm
78 pixelsize_hm = 1.18 ; % pixel size highmag tomogram in nm

```

3 Correlation from LightMicroscopy to LowMag EM image

Correlation from the original fluorescence image to an appropriate EM image containing indentifiable fiducial markers is performed using the script `martin_correlate`.

3.1 Executing the script

To execute the script and start the correlation simply run

```
martin_correlate
```

with these optional parameters in brackets

```
(fmf,emf,imf,outfileroot,fluorsel,omf,omfluor)
```

in the MATLAB command line.

```
function martin_correlate(varargin)
2
% % version MartinSchorb 130312
4 % % Copyright EMBL 2013, All rights reserved
% %
6 % =====
% % DO NOT MODIFY !!!!! use init script to set up parameters!!!
8 % =====
% %
10 % % usage is martin_correlate('beadimage','emimage','spotimage','outputfileroot','fluorsel','
% % otherimage (optional)','omfluor (optional)')
% % minimal alternative: - no parameters
12 % % - martin_correlate(outfileroot)
% % designed for correlating light and em images using fluorescent electron
14 % % dense fiducials.
% %
16 % % looks for previously picked fiducial coordinates
% %
18 % % calls cpselect for control point registration and uses cp2tform
% %
20 % % corrects for image shift between channels using bleed-thru beads
% %
22 % % uses martin_tfm_beads to suggest optimal transformation according to
% % lowest error in predictions of a single "blind" bead
24 % %
% % calls martin_corr_gui4 to check and select transformation
26 % % output overlayed transformations and predictions to screen to select
% % best transformation for the region of interest or to modify fiducials
28 % %
% % transforms fluorescence images according to the selected transformation
30 % %
% % outputs files in tif format containing (transformed) images
32 % % outputs files in tif format representing positions of picked fiducials
% % and predictions of the selected transform
34 % % (output files easily overlayed in eg imagej)
```

The following input parameters can be provided in as optional command line parameters:

1. `fmf` – path to FM image file containing fiducial information (8 or 16bit TIFF-file)
2. `emf` – path to EM image file containing visible fiducials (8 or 16bit TIFF-file)
3. `imf` – path to FM image file containing point of interest in first channel
4. `outfileroot` – directory and name base for generating output files.
5. `fluorsel` – suffix to specify type of fluorophore used ('GFP', 'YFP', ...)
6. `omf` (even more optional) – path to FM image file containing point of interest in a second channel
7. `omfluor` (even more optional) – suffix to specify second type of fluorophore used

3.2 Output and generated files

The following files are generated by the correlation script during runtime. The name base is referred to as BASE. An appended XFP refers to the fluorescent channel chosen for correlation. The selected correlation that was used is denoted by either the transform number or **all** in case the transformation based on all beads was chosen (#).

- **BASE_picked1.txt** – Plain text file containing the coordinates of fiducial pairs after subpixel fitting.

```

235         save([outfileroot,file,'.pickspots1.mat'], 'ip','bp','emf','fmf','imf','omf','slices','
            fluorsel','omfluor');
237
239     file_1 = fopen([outfileroot,file,'_picked1.txt'],'w');
    fprintf(file_1,['Picked pixel values of corresponding fluorospheres\n\nEl. Tomogram:',
        emf,'_slice:',num2str(slices.em),'_Fluorospheres:',...
        fmf,'_slice:',num2str(slices.fm),'_fluorescence image of interest:',imf,'_slice:',
            num2str(slices.im),'_other image',omf,'_slice:',num2str(slices.om),...

```

- **BASE.pickspots1.mat** – Fiducial pair coordinates, input parameters, selected fluorescent channel and clicked fluorescence spot – MATLAB format.

```

243     fprintf(file_1,'%4.2f,%4.2f_ _ _ %4.2f, _ %4.2f\n',output');
    fclose(file_1);

```

- **BASE_XFP_fluoshift.shiftcoos.mat** – Coordinates from bleed-through fiducials to determine shift in between the acquisition of images.(within **martin_chromaticshift_drift2** sub-script)

```

41         load([pathname,filename]);

```

- **BASE_XFP_#_pred.tif** – Overlay image showing the predicted and actual positions of the fiducials in EM coordinates. (size of EM-image, 16bit tiff-file)
- **BASE_XFP_#_prediction.tif** – Circle(s) marking the position(s) of the transformed coordinate(s) of the spot(s) of interest. (all images with same properties)
- **BASE_XFP_#_pred_overlay.tif** – Overlay of the prediction circle(s) and EM image
- **BASE_XFP_#_fm.tif** – Transformed fluorescent fiducial image
- **BASE_XFP_#_em.tif** – electron microscopy image
- **BASE_XFP_#_im.tif** – Transformed image of first fluorescence channel
- **BASE_XFP_#_ZFP_om.tif** (optional) – Transformed image of second fluorescence channel ZFP denotes the chosen fluorophore.
- **BASE_XFP_#_tfmed.tif** – Transformed fluorescence fiducial coordinates
- **BASE_XFP_#_pickedem.tif** – Picked EM coordinates
- **BASE_XFP_#.appltfm.mat** – MATLAB data file storing transformed coordinates, selected fluorescence channel and paths to source files for further processing
- **BASE_XFP_#_transform.log** – Plain text log file containing the source files used for correlation, the transformed spot coordinates and various information about the used transformation

EM/FM Correlation - Initial Configuration

output file name prefix

EM image file

Fiducial image file

Fluorescence image of interest

additional Fluorescence image (optional)

EM pixel size [nm]

Transformation type to use
☐ linear conformal (default)
☒ affine
☐ projective

Minimum number of beads

FM images flipped:
☒ no
☐ yes

Fluorophore:
☒ GFP
☐ RFP
☐ other

Slice number:

Slice number:

Slice number:

Slice number:

Figure 1: Initialization dialog

3.3 User interaction and key procedures

3.3.1 Initialization dialog

The first step in the correlation procedure presents a graphical interface (Fig. 1) that lets you select the input files and adjust a few parameters (all others can be set by the init script – see 2.3). The slice numbers refer to a multi-image stack, for example when providing multi-colour light microscopy images or a 3D tomogram. If provided, it will automatically use the command line parameters as values.

3.3.2 Contrast adjustment - optional

If selected in the initialization script, the contrast of each fluorescent image can be adjusted. This is done using the script `martin_contrast(image)`. The contrast can be adjusted either manually or automatically.

A file selection dialog will ask to open already existing fiducial coordinate files. If the selected images have never been used for correlating before, just close that window.

3.3.3 Rotation of the fluorescence fiducial image

If no previous fiducial coordinates could be found you can rotate the fluorescence fiducial image in 90 degree steps to simplify the initial detection of fiducial pairs. (Fig. 3) Just rotate the image until its

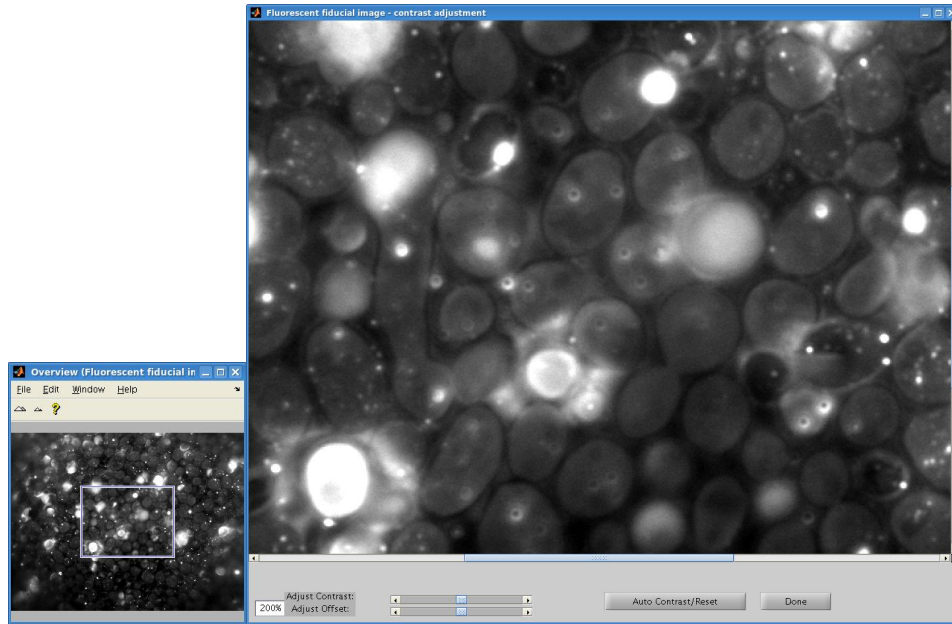


Figure 2: Contrast adjustment tool

orientation fits best the EM-image presented in the small window.

If the proper orientation is selected, click “Done”. The initial manual assignment of fiducial pairs now is done with the selected orientation. For the further processing and coordinate checks the original orientation will be used.

3.3.4 Fiducial selection

Fiducial pairs are selected in both LM and EM image using the `cpselect` tool. When an already existing coordinate file is opened, these are displayed. (Fig. 4) To continue, close the window.

3.3.5 Fiducial sub-pixel fitting and display

Fiducial positions in the light microscopy image are fitted with sub-pixel accuracy using a 2D Gaussian fit `martin_2dgaussfit`. You can specify in the init file that you want to check the fit and adjust each fiducial pair individually.

If interactive fitting is desired, a window showing the fit (Fig. 5) is presented for each of the beads. The fitted 2D-Gaussian is shown as black mesh on top of the pixel intensity visualized as a surface plot. The initially clicked position is indicated by the magenta vertical line. The following options are possible:

- ‘keep fit’ – use the fitted coordinates
- ‘keep original’ – use the original clicked coordinates without fit
- ‘select manually’ – a `cpselect` dialog showing the detailed region opens and allows you to further specify the location of the spot. Change the position of the point 1 in one of the windows (the one where the contrast fits best) and close the dialog. Then the identical fitting procedure is run on a smaller cropped region of the image around the chosen coordinates.
- ‘redo fit’ – the fit is done again using a different preset for the subtraction of linear background
- ‘kick out points’ – the current fiducial pair is removed from the dataset

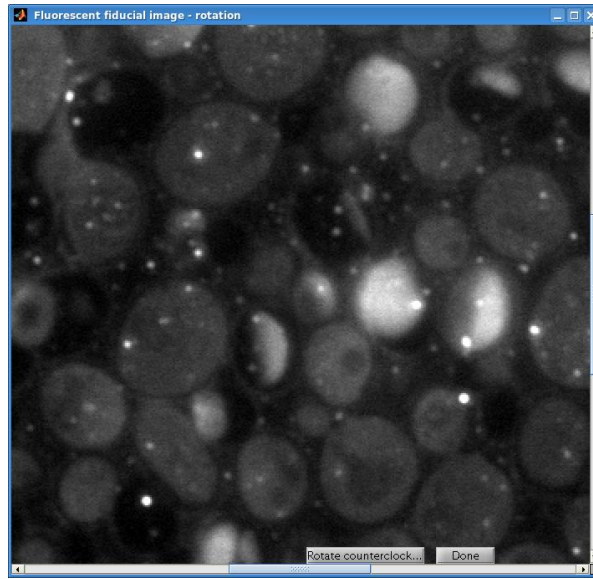


Figure 3: Image rotation window for initial correlation

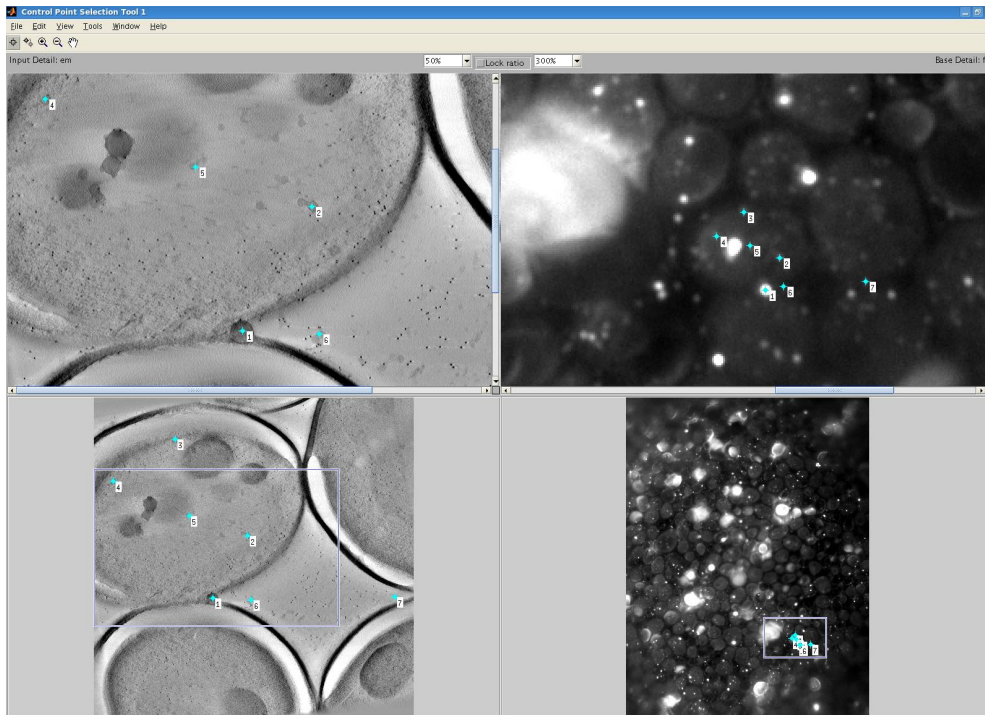


Figure 4: cpselect – graphical interface for picking and checking fiducial positions

The fitted positions are presented again using a cpselect dialog. (Fig. 4) To continue, close the window.

```

167     sixf=double(fm(floor(bp2(ispot,2))-imsir:floor(bp2(ispot,2))+imsir , floor(bp2(ispot,1))-imsir:
        floor(bp2(ispot,1))+imsir));
180     [mu,sig,Amp,check] = martin_2dgaussfit(sixf,1,fit_interactive);

        bp2(ispot,:) = floor(bp2(ispot,:)) + mu(1:2) - [1 1] - [imsir imsir];

```

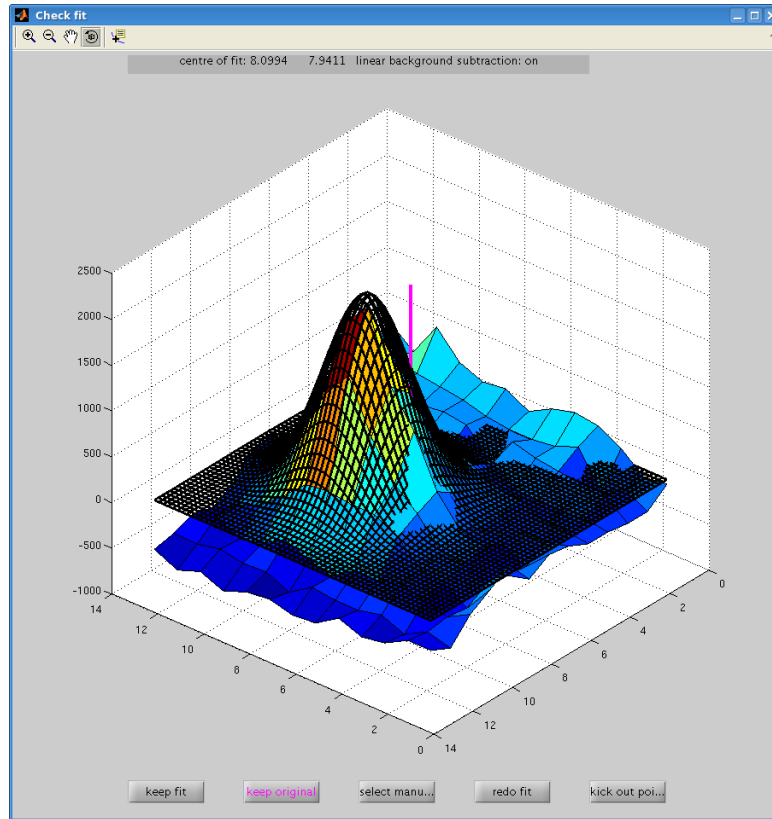


Figure 5: Interactive check of the sub-pixel fitting

3.3.6 Picking the fluorescent spot(s) of interest

In the following `cpselect` dialog, the selected fluorescence image is shown on the right, the positions indicating the fiducial markers. Click once in the right image to determine the position of the spot of interest AND once in the left image just anywhere. This click in the left image will have no effect on the correlation. `cpselect` otherwise would just not export the clicked coordinates. (Fig. 6) To continue, close the window. In case you forget to click in the left image, a reminder will be shown and you have the chance to click again.

When choosing the `multispot` option you have to select a pair of coordinates for each of the spots of interest. You can double check the intended number of spots with the clicks captured by the program.

3.3.7 Sub-pixel fitting of the fluorescent spot(s) of interest

The coordinates of the fluorescent spot(s) of interest are also determined by a gaussian fit of the image. (For details see 3.3.5, Fig. 5) The resulting coordinate(s) will be presented in a new `cpselect` window. (Fig. 7)

To continue, close the window.

```

294     sigx=double(im(floor(bpint(ispot,2))-imsir:floor(bpint(ispot,2))+imsir , floor(bpint(ispot
      ,1))-imsir:floor(bpint(ispot,1))+imsir));
296     [mu,sig,Amp,check] = martin_2dgaussfit(sigx,1,fit_interactive);
      if isnan(mu)
301         bpint1(ispot,:)=floor(bpint(ispot,:))+mu(1:2)-[1 1]-[imsir imsir];
      end

```

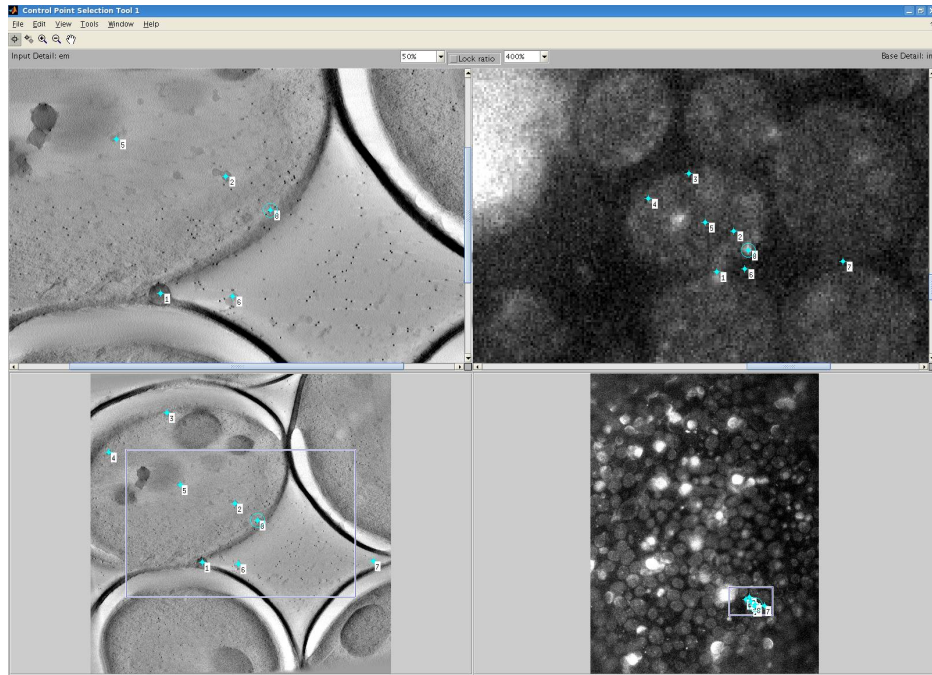


Figure 6: Selection of the fluorescent spot of interest – marked spots with nubers #8. The selected fluorescence channel image is shown on the right, the left click can be arbitrary in position.

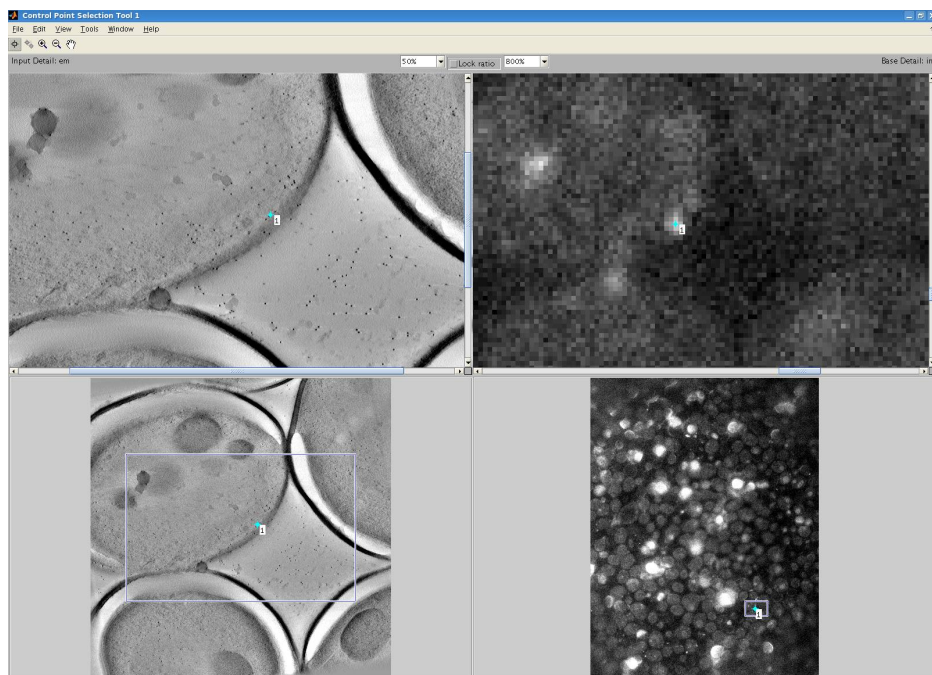


Figure 7: Positioning of the spot(s) of interest after sub-pixel fitting.

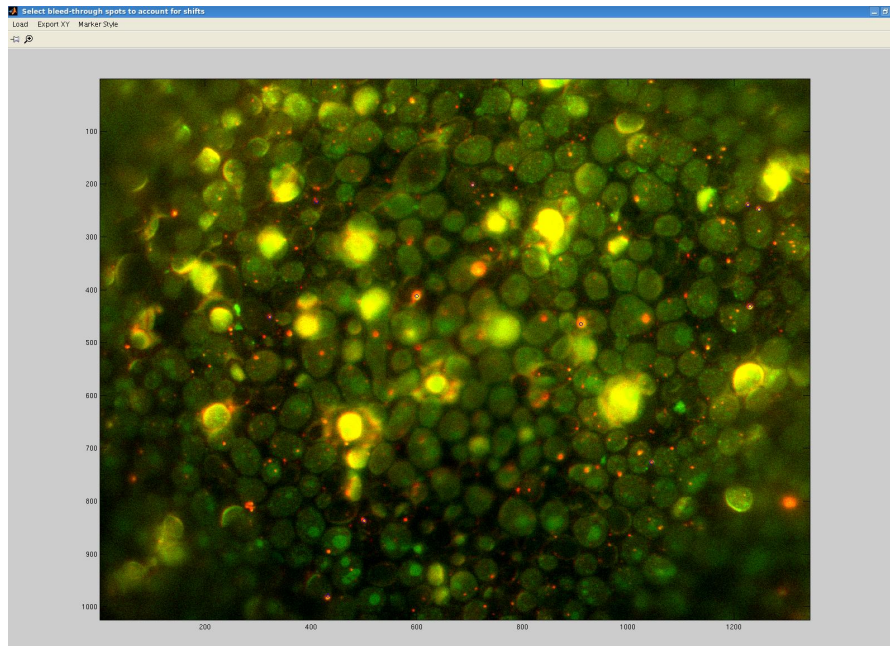


Figure 8: The tool for selecting the fiducials bleeding through to the fluorescence channel of interest. Clicked positions are indicated using blue circles.

3.3.8 Correction of stage drift in between imaging of the different fluorescence channels

The possible drift of the light microscope stage in between acquisition of the different fluorescence images is accounted for using the fiducial signal that bleeds through into the longer wavelength channels. This step can be skipped if not needed by. Just close the selection window and select “Export and Close”. To continue then press any key while having the main MATLAB window active. If you constantly want to skip this correction you can modify your configuration script accordingly.

A file selection window will ask for already existing files storing these coordinates.

In case no pre-existing shift coordinates are selected a new window will open showing an RGB-overlay of the fiducial fluorescence image in red and the selected fluorescence channel containing the signal of interest in green. (Fig. 8) Click the obvious bleeding fiducials (both channels show a spot-like signal in close proximity) and close the window, saving the positions. The entire shift correction procedure is done by:

```
339 [bluespot,fluospot]=martin_chromaticshift_drift2(fm,fm_view,im,im_view,fmboxsize,imboxsize,
    fluorsel,loc_shiftcoos,outfileroot,fit_interactive);
```

If you cannot find any obvious points, just continue without clicking positions. The shift correction will then be skipped.

To continue press any key while having the main MATLAB window active.

3.3.9 Calculation of the optimal transform and predicting fiducial and fluorescent spot EM-coordinates.

Based on the set of clicked fiducials, the most accurate linear transformation (consisting of scaling, a translation and a rotation) of coordinates is calculated using the script `martin_tfm_beads`.

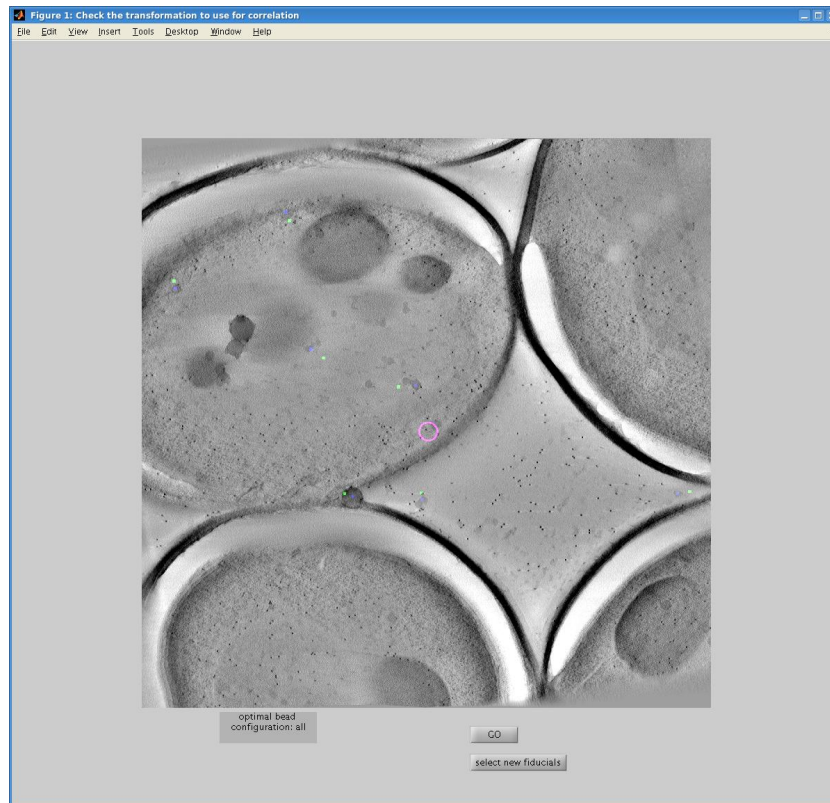


Figure 9: Preview of the coordinate transform. Fiducial selection can be modified by clicking the lower button.

```
407 % 407
[output,pickedem]=martin_tfm_beads(ip4,bp4,ipint,bpint,em,3,accuracy,init.trafo,outfileroot);
```

The resulting best coordinate transformation is applied to all coordinates and presented in a graphical output window. (Fig. 9) The beads chosen for the transformation are displayed in the bottom left. Blue positions show the clicked coordinates in the EM-image, green dots represent the transformed position of corresponding fiducial signal from the fluorescence image. To have an idea where the spot of interest is predicted, it is shown by a magenta circle. By clicking the lower button, obviously misspelled beads can be corrected, jumping to the selection step of the algorithm. (3.3.4) When clicking “GO”, all output files will be created using the presented coordinate transformation.

An overlay image showing the predicted spot(s) within the EM image as a white circle (Fig. 10) indicates a successful run of the correlating script.

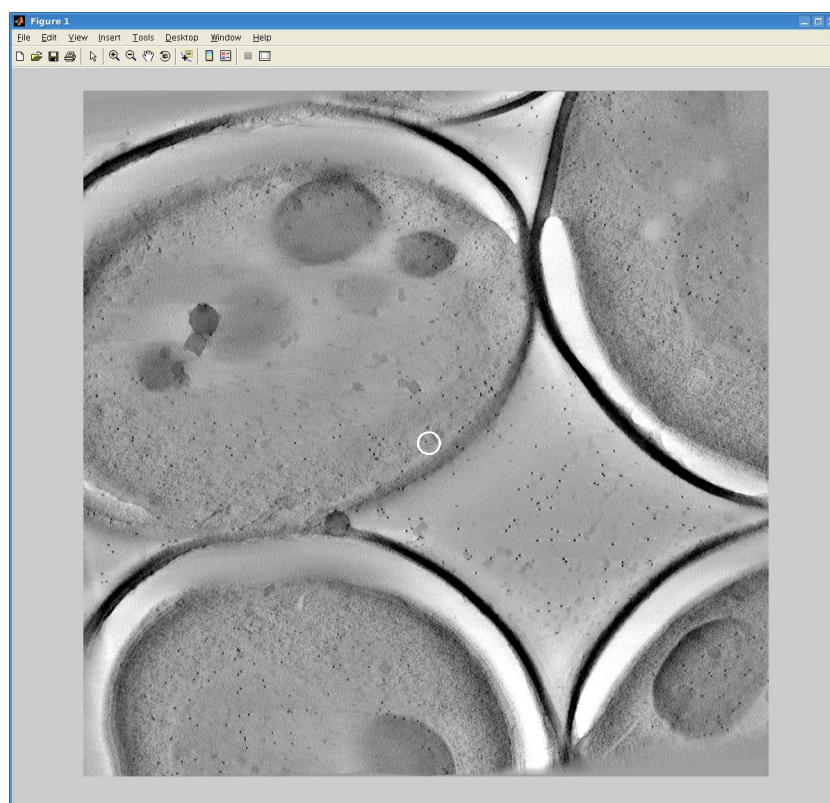


Figure 10: Performed correlation. Output position indicated by the white circle overlay.

4 Correlation from low magnification tomogram to high magnification EM image

Correlation from a low resolution EM tomogram slice to a high magnification image of the same sample containing identifiable fiducial markers such as the gold beads used for tomogram reconstruction is performed using the script `martin_LMtoHM.m`.

4.1 Executing the script

To execute the script and start the correlation simply run

```
martin_LMtoHM(hmf,smf,outfileroot)
```

in the MATLAB command line.

```
function martin_LMtoHM(hmf,smf,outfileroot,fit)
2
% % version MartinSchorb 130312
4 % % Copyright EMBL 2013, All rights reserved
%
6 %usage is martin_LMtoHM('highmaggold_image','highmag slice of interest', 'outputfileroot')
%
8 %designed for correlating em images using electron
%dense fiducials such as the gold beads used for tomogram reconstruction.
10 %
%calls cpselect for control point registration and uses cp2tform to
12 %calculate the transform
%
14 %outputs overlay registered image to screen
%outputs files in tif format containing (transformed) images
16 %outputs files in tif format representing positions of picked fiducials
%(output files easily overlayed in eg imagej)
```

It requires the following input parameters:

1. `hmf` – path to high mag. electron tomogram slice containing fiducial information (usually 2048×2048 pixel, 8 or 16bit tiff-file)
2. `smf` – path to high mag. electron tomogram slice containing a probable feature of interest that is to be checked for correlation with the fluorescence signal (usually 2048×2048 pixel, 8 or 16bit tiff-file)
3. `outfileroot` – directory and name base for generating output files.

4.2 Output and generated files

The following files are generated by the correlation script during runtime. Abbreviations of file name prefixes as in 3.2. In case you require also the transformed FM images for overlays you can activate generation of those in the config script.

- `BASE_XFP_#_XFP.lmhmcoos.mat` – Coordinates of fiducial marker pairs in both EM images
`load([pathname,filename]);`
- `BASE_XFP_#_hm.tif` – high mag. electron tomogram slice containing fiducial information
- `BASE_XFP_#_sm.tif` – high mag. electron tomogram slice containing the probable feature of interest
- `BASE_XFP_#_XFP_hm_prd_overlay.tif` – Overlay of the prediction circle and high mag EM image containing the probable feature.
- `BASE_XFP_#_XFP_hm_transform.log` – Plain text log file containing the source files used for correlation, the transformed spot coordinates and information about the used transformation



Figure 11: `cpselect` – graphical interface for picking and checking fiducial positions

4.3 User interaction and key procedures

4.3.1 Fiducial selection

A file selection dialog will ask to open the data file created by the initial correlation of the fluorescence signal. (`BASE_XFP_#.appltfm.mat`) Here you have to provide a file, otherwise the script obviously cannot transform the coordinates. A file selection dialog will ask to open already existing fiducial coordinate files for the current lowMag to highMag correlation. If the selected images have never been used for correlating before, just close that window.

Fiducial pairs are selected in both LM and HM image using the `cpselect` tool. When an already existing coordinate file is opened, these are displayed. (Fig. 11) To continue, close the window.

4.3.2 Correlation

Coordinates are transformed using all marked fiducials and the output files are written. An overlay image showing the predicted spot within the highMag image containing the feature of interest (Fig. 12) indicates a successful run of the correlating script.



Figure 12: Performed correlation. Output position indicated by the white circle overlay.

5 Useful hints and tricks

5.1 Overlay of correlated images

The correlation procedure writes a variety of images in tiff format that have all be aligned to each other according to the transformation used for correlation. (see 3.2) These include transformed fluorescence images as well as prediction circles or fiducial positions. You can easily superimpose these images to the EM image or even the tomogram.

ImageJ (version>1.41)

Open both the EM image and the image you want to superimpose in ImageJ/Fiji. In Image>Color>Merge Channels select the colour channel of choice for the overlay(s) and grey for the EM image. You can then adjust the contrast for each channel by selecting the colour using the mouse wheel or the slider at the bottom of the image and the contrast adjustment tool.

Amira

Open your original tomogram (beware of the file format issue in Amira) and create an OrthoSlice. Make sure the perspective is set to parallel. This can be done by clicking the tool button showing an eye and either parallel or diverging beams. 

Open the image you want to superimpose, accept the parameters suggested (voxel size should be 1,1,1) and create an OrthoSlice. Select the transparency to alpha. If you want you can also adjust the color mapping.